

Aesthi Peel™ — Aesthetic Microdermabrasion Unit

Easy to use

Sterile disposable consumables

Patented 45°, v- shaped, light-weight handpiece

Portable and light weight

Power adjustment features

Suitable for:

Facial and body peeling

Aging skin

Dull, congested and dry skin


Enlarged pores

Superficial rhytids

Minor scars

Localized hyperpigmentation

Pre-treatment for further applications (Laser and Chemical Peel)

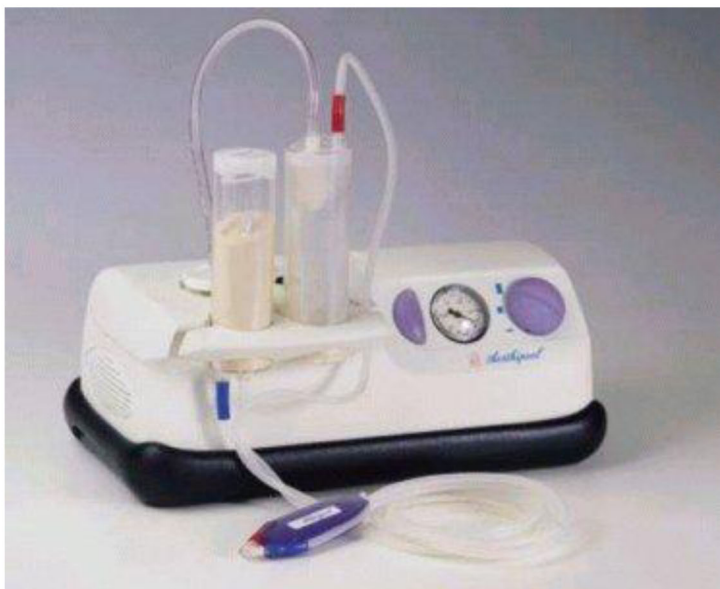

The Aesthi Peel™ microdermabrasion unit is a pure aesthetic device manufactured by Mattioli Engineering.

Aesthi Peel™ performs a progressive and controlled exfoliation of the skin.

Aesthi Peel™ uses the Mattioli Engineering patented 45 degree, V-shaped handpiece and other patented, disposable supplies. It is transportable and clog-resistant.

This specially engineered handpiece allows crystals to flow in a uniform and effective way, abrading the skin in a comfortable fashion.

Aesthi Peel™ allows the users to work with the sterilized disposable delivery kit (i.e. bottles, corundum, tubing kit, and handpiece) in order to perform a clean and safe operation. Cross contamination risks are completely avoided.

The light, easy-to-use accessories, as well as the preset levels of vacuum power, make the Aesthi Peel the ideal solution for superficial and medium peeling procedures.

ACCESSORIES AND DISPOSABLES



*Pre-Packed, Pre-Sterilized,
disposable delivery Kit,
accessories and corundum
30000385 Corundum
bottle set
30010010 Disposable
Tubing kit
30010300 Air Filters*



*45 Deg. V-Shaped
Patented disposable
abrading tip
10010142*



*Carrying case w/trolley
available upon request
CARRYING CASE*

TECHNICAL SPECIFICATIONS

POWER SUPPLY:

Supply Voltage: 110 V, single phase
Frequency: 60 Hz
Power absorption: 200 VA
Fuses: 2 x 1 AT at 110V, 60Hz
Class (according to IEC 601-1) 1 type BF

PERFORMANCES:

Suction power: from 0.2 to 0.7 bar
(below 10000 mt altitude)

Abrading media:

High purity medical grade corundum powder

PHYSICAL:

Dimension : 440x300x175 mm
173x118x69 mm
Weight 8 kg (17,6lb)

STERILIZATION


Sterilization is of paramount importance when the operator is dealing with different patients and skin conditions.

Cross contamination and debris may affect the patient if the microdermabrasion device is not perfectly cleaned and sterile.

Only Mattioli Engineering guarantees safe and sterile operations through its disposable accessories and consumables.

Mattioli Engineering - A new science in Skin Caretm

Mattioli Engineering Italia SPA
VIALE MACHIAVELLI 2/A
50125 FIRENZE—ITALY
Phone.: +39 055 221603
Fax: +39 055 221735
Mattioli@ftbcc.it

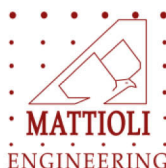

Saratoga Aesthetics
6909 Sandy Shore Road
Columbia, SC 29206
1.888.785.7546 phone
1.803.788.7776 fax
info@saratogaaesthetics.com
www.SaratogaAesthetics.com

DISTRIBUTED BY: SARATOGA AESTHETICS 1-888-785-7546